

TEST REPORT

扫码下载报告



VERIFICATION WEBSITE: www.gttc.net.cn

VERIFICATION CODE: MHHC-0461-44

ISSUE DATE: 2020-05-21



No: 200100001

APPLICANT: HUIZHOU YIMIAOAN LABOR INSURANCE SUPPLIES CO., LTD
ADDRESS: THE 3RD FLOOR, BUILDING NO. F, XIawei VILLAGE FACTORY (HUIZHOU LEISHI INDUSTRIAL PARK), RUHU TOWN, HUICHENG DISTRICT, HUIZHOU CITY

FORMATION CONFIRMED BY APPLICANT:

95 MASK

QUANTITY: THIRTY PIECES

DATE RECEIVED/DATE TEST STARTED: 2020-04-29

CONCLUSION:

FILTRATION EFFICIENCY TO NaCl PARTICULATE MATTER	M
RESPIRATORY RESISTANCE	M
RESPIRATORY RESISTANCE	M
PERMEABILITY	M
HEAD BAND	M

NOTE: "M" -MEET THE STANDARD'S REQUIREMENT "F" -FAIL TO MEET THE STANDARD'S REQUIREMENT
"---" -NO COMMENT

REMARK:

THIS REPORT IS THE ENGLISH TRANSLATION VERSION OF THE REPORT 200095008.
ALL THE TESTED ITEMS ARE TESTED UNDER THE STANDARD CONDITION (EXCEPT FOR INDICATION).
ALL COPIES OF THE REPORT ARE VALID ONLY RE-STAMPED.
THE EXPERIMENT WAS CARRIED OUT AT No. 1, ZHUJIANG ROAD, PANYU DISTRICT, GUANGZHOU, GUANGDONG, P. R. CHINA.


APPROVED BY:
Ma ENGINEER

马楠



TEST REPORT

No: 200100001

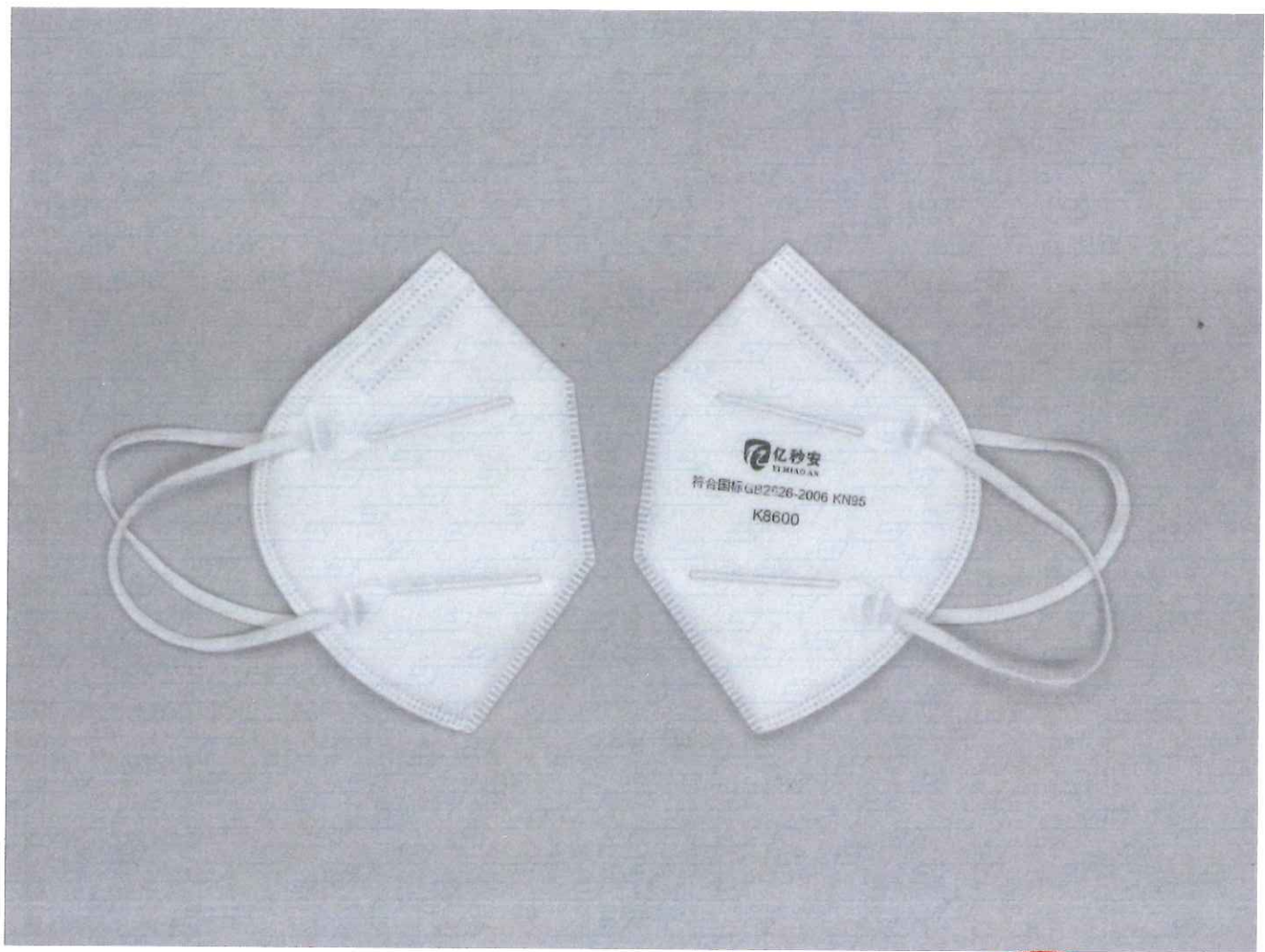

TEST REPORT

No: 200100001

FILTRATION EFFICIENCY TO NaCl PARTICULATE MATTER (%)

(GB 2626-2006 6.3, AIR FLOW: 85L/min, AEROSOL: NaCl, AEROSOL CONCENTRATION: 15mg/m³, TEMP: 22.5°C RH: 36.2%)

FILTRATION EFFICIENCY:

UNTREATED SAMPLE

1# 99.407
2# 99.356
3# 99.725
4# 99.418
5# 99.329
6# 99.146
7# 99.627
8# 99.270
9# 99.535
10# 99.203

CONDITIONING TREATED

1# 99.078
2# 99.720
3# 99.468
4# 99.507
5# 99.276

REQUIREMENT

FILTRATION EFFICIENCY:

≥95.0
(KN95)
(GB 2626-2006)

INSPIRATORY RESISTANCE (Pa)

(GB 2626-2006 6.5, HEAD SIZE: MEDIUM)

UNTREATED SAMPLE:

1# 165.3
2# 151.4

PRETREATMENT SAMPLE:

1# 153.5
2# 143.5

REQUIREMENT

≤350
(GB 2626-2006)

EXPIRATORY RESISTANCE (Pa)

(GB 2626-2006 6.6, HEAD SIZE: MEDIUM)

UNTREATED SAMPLE:

1# 152.0
2# 141.4

PRETREATMENT SAMPLE:

1# 138.6
2# 132.1

REQUIREMENT

≤250
(GB 2626-2006)

FLAMMABILITY (s)

(GB 2626-2006 6.15)

AFTERFLAME TIME

UNTREATED SAMPLE

1# 0.0
2# 0.0

CONDITIONING TREATED

3# 0.0
4# 0.0

REQUIREMENT

AFTERFLAME TIME
≤5
(GB 2626-2006)



TEST REPORT

No:200100001

HEAD BAND

(GB 2626-2006 6.11)

UNTREATED SAMPLE:

1# PASS

CONDITIONING TREATED:

1# PASS

REQUIREMENT

EACH HEAD BAND, BUCKLE AND OTHER ADJUST
PARTS OF MASK SHOULD NOT SLIP OR BREAK

UNDER 10N TENSION FOR 10S.

(GB 2626-2006)



——End of Report——

PAGE 4 OF 4