


**TECHNICAL DATA SHEET**  
**Self-Seal Sterilization Pouch**

**Intended Use**

The self-sealing sterilization pouches are intended to provide health care workers with an effective method to enclose devices intended for sterilization in steam autoclaves or via Ethylene Oxide (EO). The pouch's process indicators are designed to indicate to the user that the pouch has undergone either a Steam or EO sterilization process. Single use.

Characteristic		Specification		
Material Composition	Film	PET and CPP laminated film		
	Paper	Medical grade kraft paper		
Description	Laminated film sealed to paper with two sets of indicators and double-sided adhesive tape.			
Process Indicator Exposure	Steam: Indicators turn pink to brown when processed for 10 min. at 121°C or 2 min. at 134°C EO: Indicators turn blue to gold when processed for 20 min at 54° C, 60% RH, or 30 min at 30° C, 60% RH and EO concentration of 600 mg/L for both options.			
Color	Refer to Appendix 1			
Dimension of the film	Size	Width (mm)	Length (mm)	
	2.25" × 4"	57	101	
	2.75" × 9"	70	225	
	3.5" × 5.25"	89	137	
	3.5" × 9"	89	225	
	5.25" × 10"	144	251	
	7.5" × 13"	191	330	
	9" × 11"	229	279	

**Product Image**

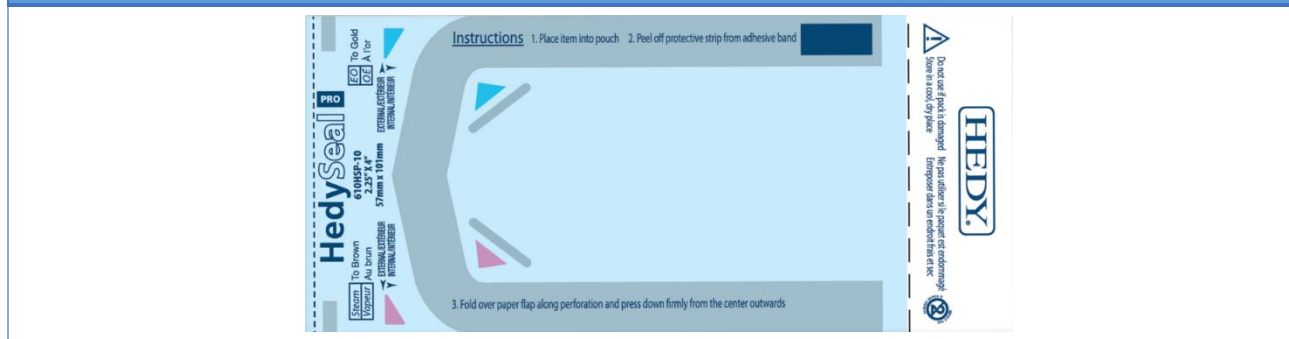



#### Additional Information

Storage Condition	Country of Origin
Keep away from light, heat, and humidity. Store at temperature between 15 – 35°C (59 – 95°F) with a relative humidity of <70%	Canada

#### Appendix 1 List of Item Codes and Attributes

## Appendix 1 (List of items)

*AMD Medicom Inc.*

Item Code	Description	Packaging Information	Size	Color	Treatment
610HSP-10	HedySeal Pro Class 4 Self-Sealing Sterilization Pouches	200 Pouches/Box, 12 Boxes/Case, 2400 Pouches/Case	2.25" x 4"	Clear Film	Non-Sterile
610HSP-20	HedySeal Pro Class 4 Self-Sealing Sterilization Pouches	200 Pouches/Box, 6 Boxes/Case, 1200 Pouches/Case	2.75" x 9"	Clear Film	Non-Sterile
610HSP-30	HedySeal Pro Class 4 Self-Sealing Sterilization Pouches	200 Pouches/Box, 6 Boxes/Case, 1200 Pouches/Case	3.5" x 5.25"	Clear Film	Non-Sterile
610HSP-40	HedySeal Pro Class 4 Self-Sealing Sterilization Pouches	200 Pouches/Box, 6 Boxes/Case, 1200 Pouches/Case	3.5" x 9"	Clear Film	Non-Sterile
610HSP-50	HedySeal Pro Class 4 Self-Sealing Sterilization Pouches	200 Pouches/Box, 6 Boxes/Case, 1200 Pouches/Case	5.25" x 10"	Clear Film	Non-Sterile
610HSP-60	HedySeal Pro Class 4 Self-Sealing Sterilization Pouches	200 Pouches/Box, 5 Boxes/Case, 1000 Pouches/Case	7.5" x 13"	Clear Film	Non-Sterile
HP-0911	Hedy Pouch-ette Self-Sealing Sterilization Pouches	200 Pouches/Box, 5 Boxes/Case, 1000 Pouches/Case	9" x 11"	Blue Film	Non-Sterile